


Product Specification Sheet

| | |
|------------------------------|--|
| Product Name | pLenti-MCS-PGK-HygR Lentiviral Vector |
| Description | <p>Lentivirus vector based on the human immunodeficiency virus-1 (HIV-1) has become a promising vector for gene transfer studies. The advantageous feature of lentivirus vector is the ability of gene transfer and integration into dividing and nondividing cells. The pseudotyped envelope with vesicular stomatitis virus envelope G (VSV-G) protein broadens the target cell range. Lentiviral vectors have been shown to deliver genes to neurons, lymphocytes and macrophages, cell types that previous retrovirus vectors could not be used. Lentiviral vectors have also proven to be effective in transducing brain, liver, muscle, and retina in vivo without toxicity or immune responses. Recently, the lentivirus system is widely used to integrate siRNA efficiently in a wide variety of cell lines and primary cells both in vitro and in vivo. Lentivirus particles are produced from 293T cells through transient transfection of plasmids that encode for the components of the virion. Our third generation lentiviral systems have been designed for increased researcher safety. pLenti-MCS-PGK-HygR Lentiviral Expression Vector contains the MCS for cloning of custom Promoter-GOI cassette and HygR gene driven by PGK promoter. The PGK-HygR cassette permits hygromycin selection.</p> |
| Catalog Number | LV003 |
| Size | 10 µg at 0.5 µg/µL in TE |
| Shipping | Room temperature |
| Storage and Stability | Store at -20°C immediately upon receipt. This product is stable for 6 months when stored as directed. |
| Quality Control | This plasmid is sequence verified. |
| Safety Precaution | Remember that you will be working with samples containing infectious virus. Follow the recommended NIH guidelines for all materials containing BSL-2 organisms. The ALSTEM Lentiviral Expression System is designed to minimize the chance of generating replication-competent lentivirus, but precautions should still be taken to avoid direct contact with viral supernatants. |
| Restricted Use | For Research Use Only. Not for use in diagnostic or therapeutic procedures. |

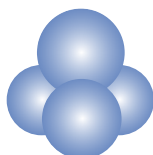

ALSTEM, INC

2600 Hilltop Drive, BLDG B, STE C328, Richmond, CA 94806

Tel: (510) 708-0096


Fax: (866) 605-8766

www.alstembio.com

info@alstembio.com

Vector Information


This is a lentiviral expression vector that contains all the elements necessary for efficient and high yield viral production. The MCS provides opportunities for cloning of custom Promoter-GOI cassette and PGK-HygR cassette allows for selection of transduced cells. This vector can be used for stable cell line generation.

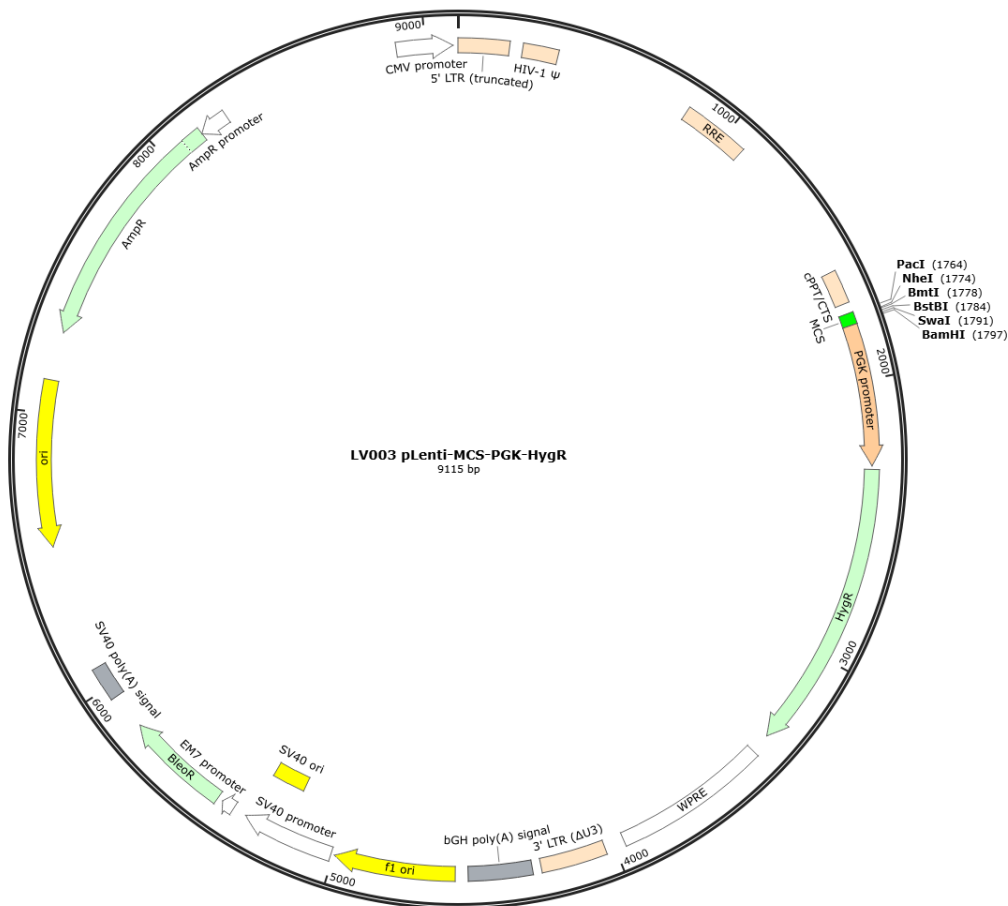

Note: Bacterial culture of pLenti vectors should be done in medium containing 50 µg/mL Ampicillin. For maximal plasmid yield and quality, we recommend Stbl3 competent cells (Invitrogen).

IMPORTANT NOTICE

Store the vial at -20°C immediately upon receipt.

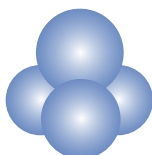

ALSTEM, INC

2600 Hilltop Drive, BLDG B, STE C328, Richmond, CA 94806

Tel: (510) 708-0096

Fax: (866) 605-8766

www.alstembio.com

info@alstembio.com